

B-SAFE™ SEATED SLING WITH HEAD SUPPORT

Available with clip or loop attachments

B-SAFE™ Seated Sling Single Patient Use slings support the whole body including the head. They are designed for use on one patient during their stay and should not be redistributed to other patients. With built-in head support and padded legs, the *B-SAFE™ Seated Sling* ensures superior patient comfort during transfers. The *B-SAFE™ Seated Sling* can be used for bed, chair and floor lifts.

- Slings are available with loop or clip attachments to fit clip or loop style lifters
- Available in weight capacities from 600lbs to 1000lbs
- Made of 70% Polyester / 30% Viscose with Polypropylene Straps
- **Always ensure you have read the Sling Application Guide before use**

BENEFITS FOR STAFF AND PATIENTS

- Can be used in pools, bathing and showers
- Soft Poly / Viscose to maintain skin integrity
- Breathable - this reduces sweating to help reduce moisture related skin integrity (see air permeable report)
- Absorbent - fluids don't run off
- Minimal stretch under load weight
- Can be used during XRAY / MRI



Identification label

The sewn-on label has space to write the patient's name to clearly identify the sling as theirs, to reduce any risk of cross-infection.

Wash indicator

The *B-SAFE™ Seated Sling* is NOT designed to be washed. If the sling is accidentally washed, the wash indicator will change labels so that the care provider is aware the sling not be used further.



Expected lifetime

This sling can be used as long as it does not display any signs of wear or damage prior to each use, however we recommend that you replace the sling after 6 months of use

Easy identification


Size is indicated by the colour of the binding

Extra strength

Unlike some disposable slings, the *B-SAFE™ Seated Sling* is Polyester Viscose for extra strength

Padded legs for comfort

Attachment types

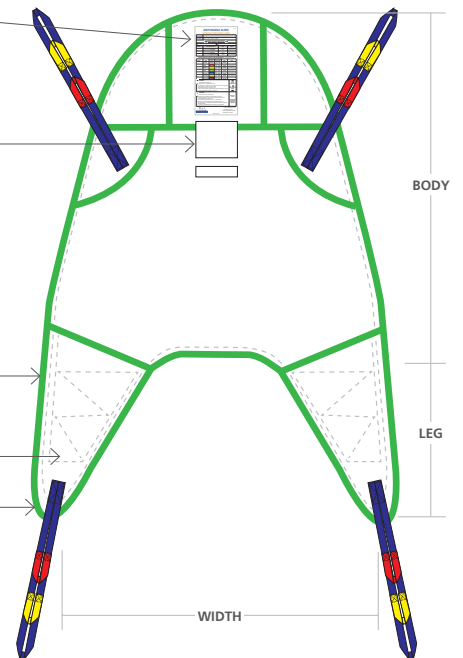

B-SAFE™ REPOSITIONING SLING SINGLE PATIENT USE

Available with clip or loop attachments

B-SAFE™ Repositioning Single Patient Use is a full bed-length sling, used with 4 point spreader bars, that allows the patient or resident to remain in a supine position for transfer or in-bed repositioning. It can also be used for turning. They are designed for use on one patient during their stay and should not be redistributed to other patients.

This sling is compatible with hook style lifters (contact *Bridge Healthcare* for further information – it is not compatible with clip style lifters).


They are most often used as a replacement for reusable slings:

- In high infection risk areas where the stitching presents a risk
- Where wash costs and lead times are high
- To reduce high initial capital expenditures where budgets don't allow for it
- During outbreaks of known or unknown infectious disease or Gastroenteritis



Reinforced edges for patient comfort

Straps are sewn and reinforced onto the edge so the straps are not under the patient and pushing into their skin.

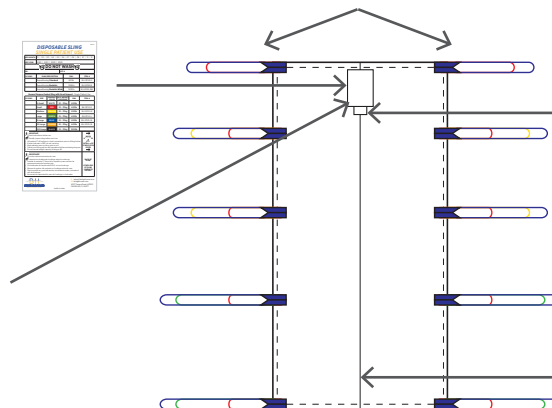

Identification label

The sewn-on label has space to write the patient's name to clearly identify the sling as theirs, to reduce any risk of cross-infection.


DO NOT WASH

The sewn on label also clearly highlights **DO NOT WASH** with the international symbol



Wash indicator

The sling is **NOT** designed to be washed. If the sling is accidentally washed or heavily soiled, the wash indicator will turn red so that the care provider is aware the sling is not to be further used.

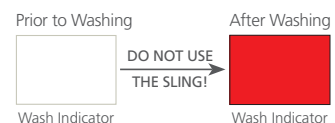

Line for patient alignment (Bariatric only)

Always ensure you have read the sling application guide before use



To request a free sample, send an email to sales@bridgehcusa.com